

H201-PRIOR-NZ500

Stage compatibility: screws onto Prior NZ250/400/500 - comes with Koehler lid

Magnetic interchangeable inserts allow to use:

- 35/60 mm Petri
- 1x3 in. / 1x2 in. chambered slides
- 6-12-24-48-96 MW plates

Features:

- Koehler lid for working under Koehler illumination
- Optional lids available for Laser Safety, Injection during imaging
- Flexible chamber height thanks to the removable chamber riser required to fit standard multiwell plates or to perform perfusion
- Magnets allow easy interchange of chamber inserts
- Magnetic locks hold Petri and Slides and MW in correct position inside the chamber




AVAILABLE INSERTS

H201-MW-HOLDER-NZ500	Multiwell plates holder
1x35-M	#1 35mm Petri-dish
1x35-PRIOR-NZ500	#1 35mm Petri-dish
1xGS-PRIOR-N-NZ500	#1 1x3in. chamberslide
1xLABTEK-PRIOR-NZ500	#1 Lab-Tek 1x2 in. chambered cover glass
1xLABTEK-II-PRIOR-NZ500	#1 Lab-Tek II 1x2 in. chambered cover glass
1x60-PRIOR-NZ500	#1 60mm Petri-dish
1xT25-PRIOR-NZ500	#1 Nunc and Greiner T25 flask
2x35-PRIOR-NZ500	#2 35mm Petri-dish
2xGS-PRIOR-NZ500	#2 1x3in. chamberslides
2xLABTEK-PRIOR-NZ500	#2 Lab-Tek 1x2 in. chambered cover glass
2xLABTEK-II-PRIOR-NZ500	#2 Lab-Tek II 1x2 in. chambered cover glass
2x60-PRIOR-NZ500	#2 60mm Petri-dish
4x35-PRIOR-NZ500	#4 35mm Petri-dish
GS35-PRIOR-NZ500	#1 1x3in. chamberslide and #2 35mm Petri-dish
LABTEK-35-PRIOR-NZ500	#1 Lab-Tek 1x2 in. chambered cover glass and #2 35mm Petri-dish
LABTEK-II-35-PRIOR-NZ500	#1 Lab-Tek II 1x2 in. chambered cover glass and #2 35mm Petri-dish
LABTEK-60-PRIOR-NZ500	#1 Lab-Tek 1x2 in. chambered cover glass and #1 60mm Petri-dish
LABTEK-II-60-PRIOR-NZ500	#1 Lab-Tek II 1x2 in. chambered cover glass and #1 60mm Petri-dish
GS60-PRIOR-NZ500	#1 1x3in. chamberslide and #1 60mm Petri-dish
6035-PRIOR-NZ500	#1 60mm Petri-dish and #1 35mm Petri-dish

AVAILABLE LIDS

H201-KOEHLER-LID-PRIOR-NZ500	Glass lid that reduces the height of the chamber to 22 mm. It allows to work under Koehler illumination.
H201-INJECTION-LID-PRIOR-NZ500	Glass lid with two small openings (sealed with flexible plastic) to allow injection or permanent access to the specimen. Note: compatible only with specimen holder 1x35
H201-LASER-INTERLOCK-LID-PRIOR-NZ500	Glass lid with a safety switch that connects to the laser box. The switch turns off the laser, when the lid is lifted Note: compatible only with specimen holder 1x35

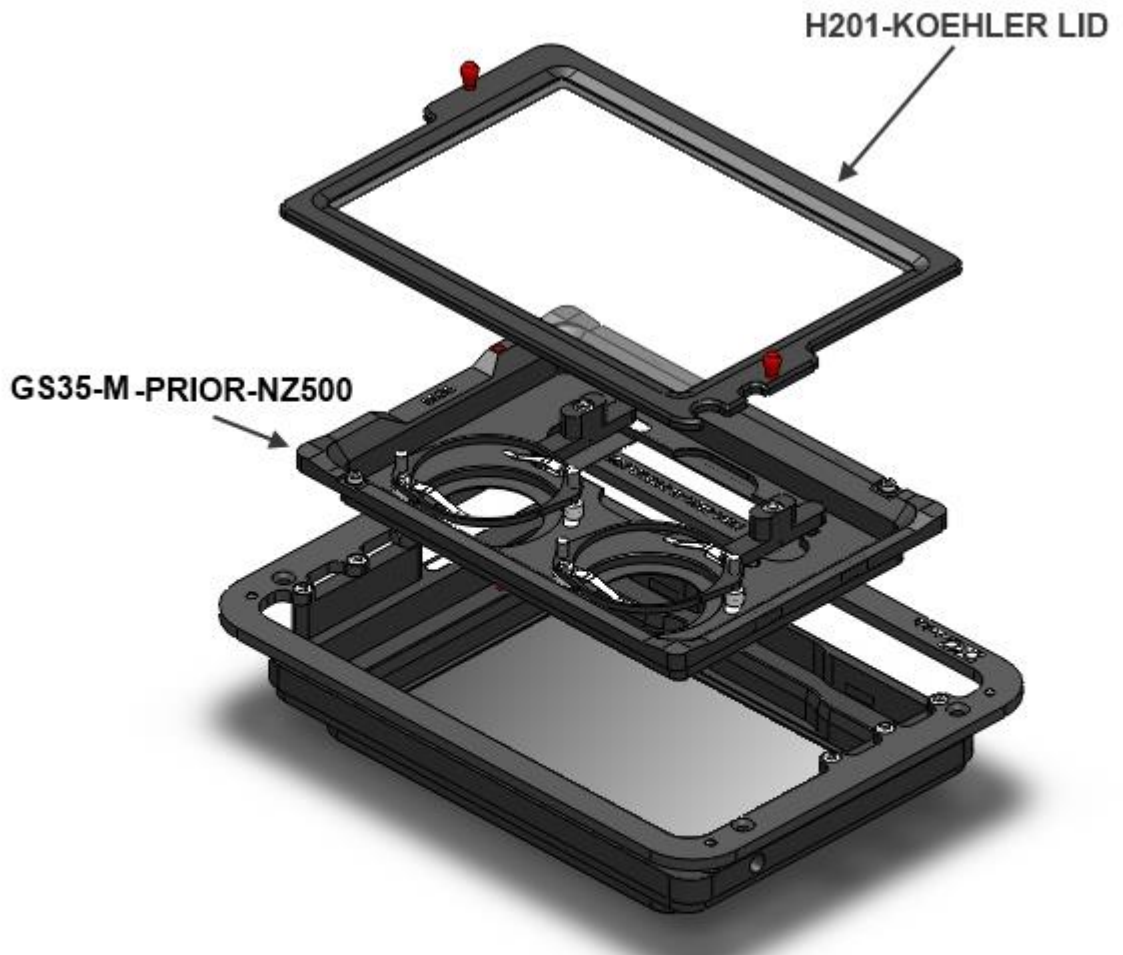

Ordering Info
H201-PRIOR-NZ500
H201-KOEHLER-LID
GS35-PRIOR-NZ500

H201-PRIOR-NZ500 with Koehler Lid



H201-KOEHLER LID-PRIOR-NZ500
Included as standard






H201-LASER INTERLOCK LID-PRIOR-NZ500



H201-INJECTION-LID-PRIOR-NZ500


H201-PRIOR-NZ500 - Available Lids

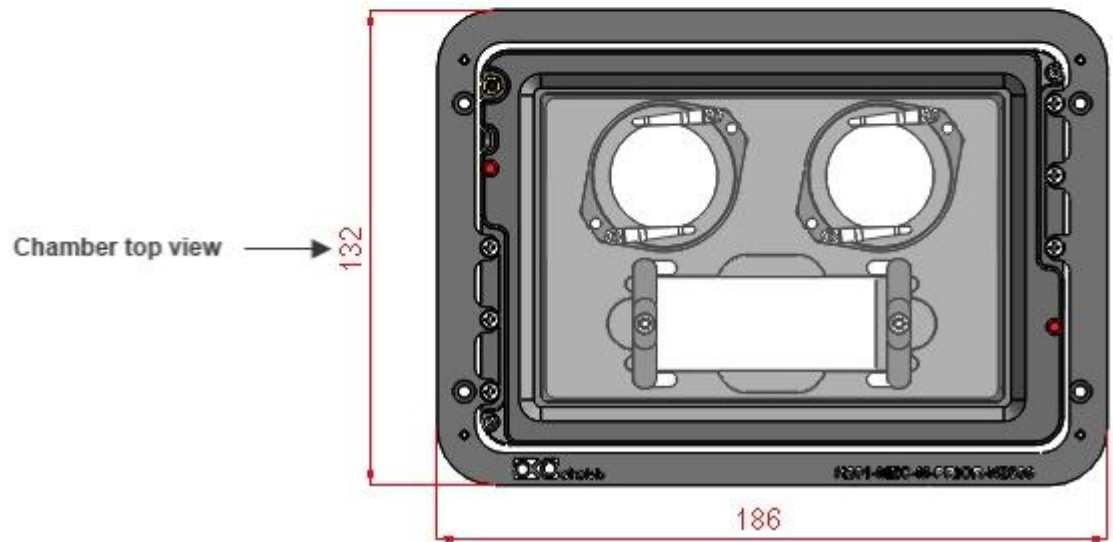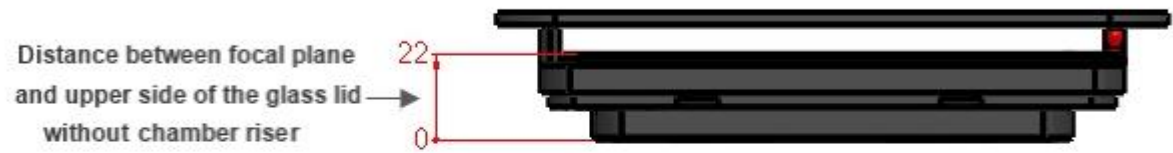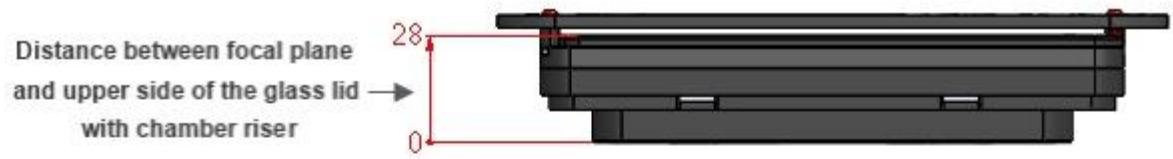


H201-PRIOR-NZ500 - Dimensions




H201-PRIOR-NZ500 - Magnetic Inserts

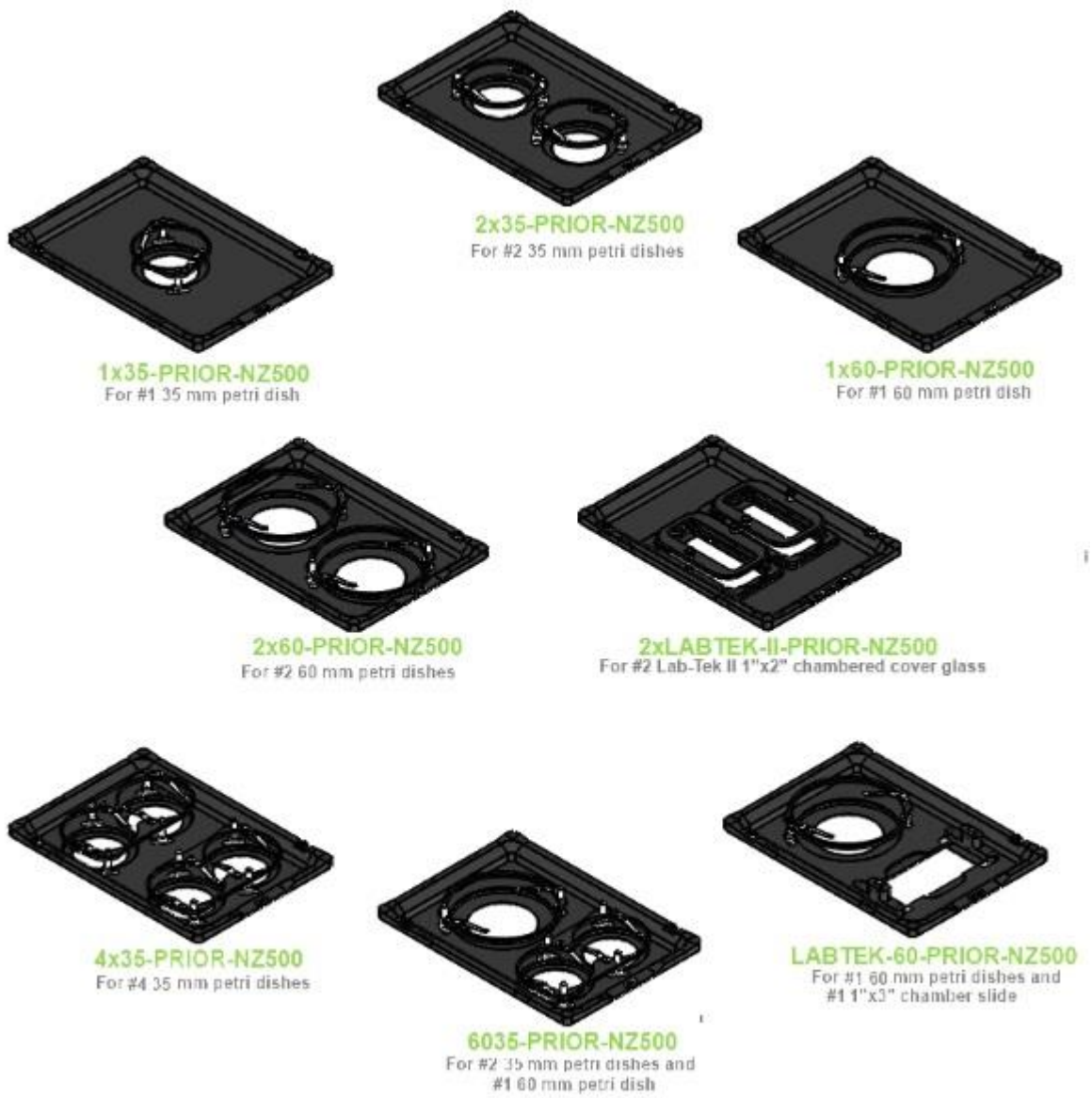

H201-PRIOR-NZ500 - Magnetic Inserts

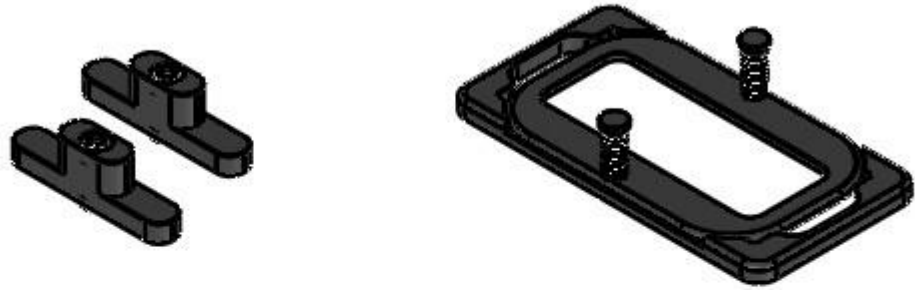

**Magnetic Locks for chamber slides and chambered cover glass.
Included in the corresponding chamber insert**

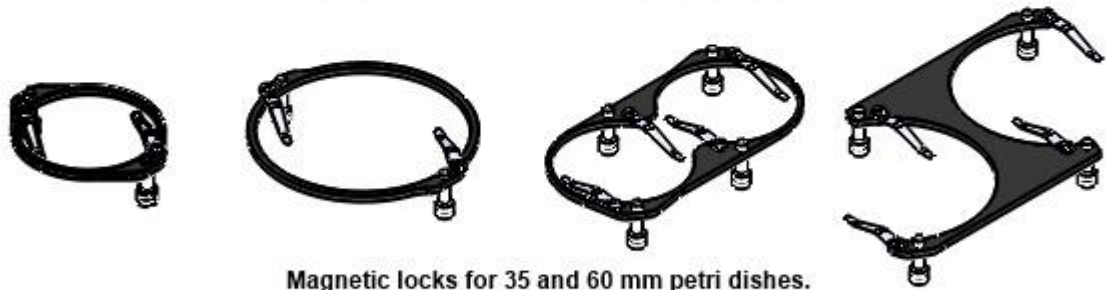

**Magnetic locks for 35 and 60 mm petri dishes.
Included in the corresponding chamber insert**

H201-PRIOR-NZ500 - Magnetic sample locks