





FOR INFORMATION ONLY
NOT FOR MANUFACTURING PURPOSES

DRAWING PROJECTION				THORLABS www.thorlabs.com	
NAME		DATE		DODT GRADIENT CONTRAST MODULE	
DRAWN	TB	07/JUN/16			
APPROVAL	PN	07/JUN/16		MATERIAL	
COPYRIGHT © 2016 BY THORLABS				N/A	
VALUES IN PARENTHESIS ARE CALCULATED AND MAY CONTAIN ROUND OFF ERRORS				ITEM #	REV
				WFA1100	A
				APPROX WEIGHT	
				1.6 kg	


BOM	
①	DIFFUSER SLIDER
②	COVER
③	M2.5 x 6mm FLAT HEAD MACHINE SCREWS

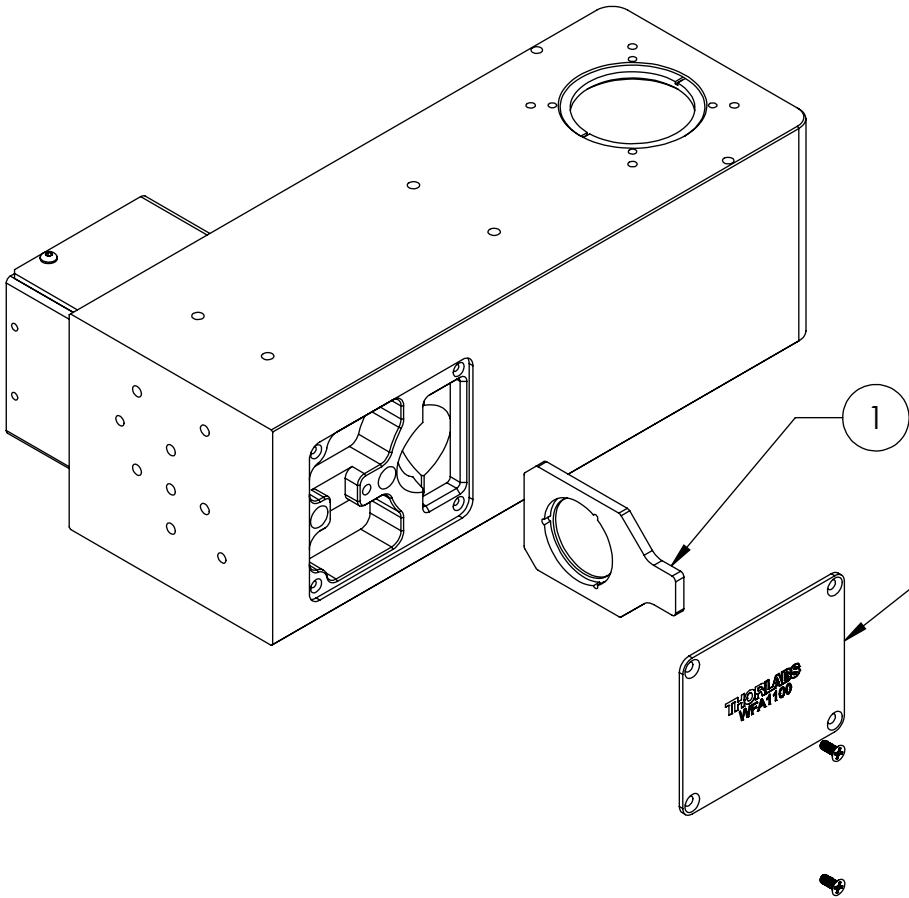

NOTE: REMOVE DIFFUSER TO RUN DODT CONTRAST.
TO RUN BRIGHTFIELD USING THIS MODULE,
KEEP DIFFUSER IN, AND REMOVE ROTATION MOUNT
SHOWN ON THE NEXT PAGE


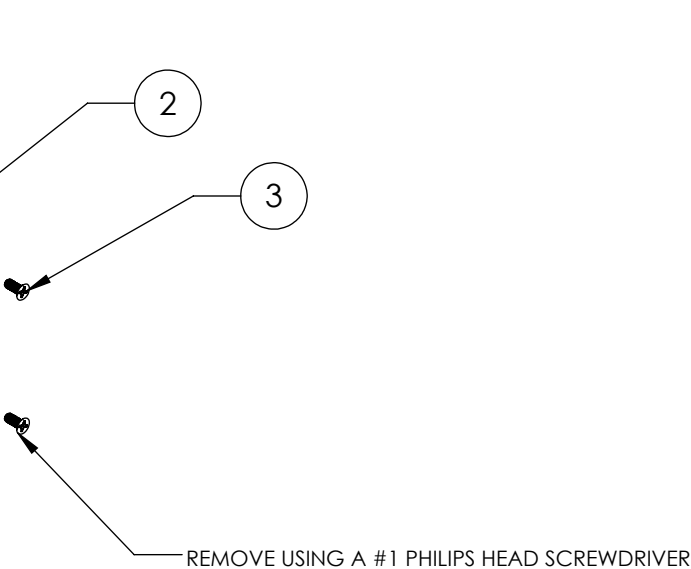


FOR INFORMATION ONLY
NOT FOR MANUFACTURING PURPOSES


DRAWING PROJECTION			 www.thorlabs.com	
	NAME	DATE		
DRAWN	TB	07/JUN/16	DODT GRADIENT CONTRAST MODULE	
APPROVAL	PN	07/JUN/16	MATERIAL	
COPYRIGHT © 2016 BY THORLABS			N/A	
VALUES IN PARENTHESIS ARE CALCULATED AND MAY CONTAIN ROUND OFF ERRORS			ITEM #	APPROX WEIGHT
			WFA1100	1.6 kg

REV
A


STEP 1:
REMOVE COVER AND
2.5mm x 6 SCREWS
USING A 1.5mm HEX

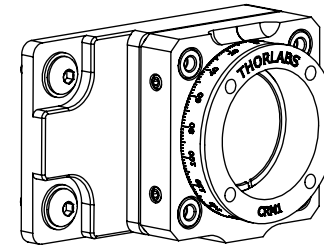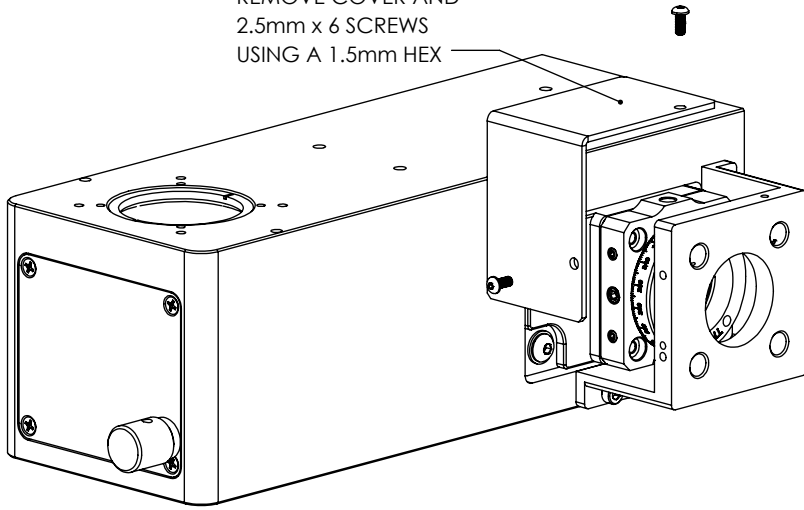

DETAIL A
SCALE 2 : 3

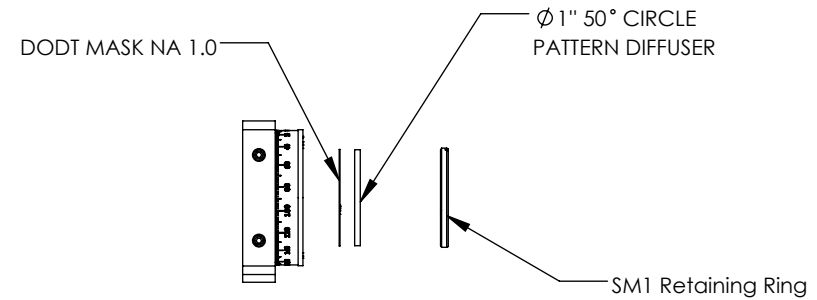


STEP 2:
REMOVE ROTATION MOUNT
TO REPLACE MASK. SEE DETAIL A

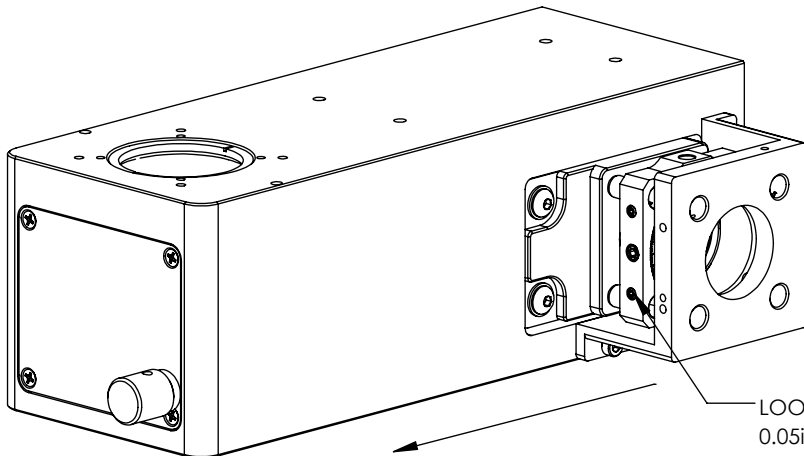

LOOSEN LOCKING SETSCREW
0.05in [1.3mm] HEX
2 PLACES (BOTH SIDES)
AND SLIDE OUT ROTATION MOUNT

THE FOLLOWING ADDITIONAL MASKS ARE INCLUDED:

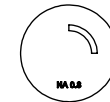


DODT MASK NA 0.8



DODT MASK NA 0.65





DODT MASK NA 0.5

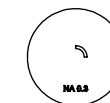

DODT MASK NA 0.3

FOR INFORMATION ONLY
NOT FOR MANUFACTURING PURPOSES

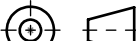DRAWING PROJECTION			 www.thorlabs.com	
NAME		DATE	MASK AND ACCESSORIES FOR DODT	
DRAWN	TB	07/JUN/16		
APPROVAL	PN	07/JUN/16	MATERIAL	
COPYRIGHT © 2016 BY THORLABS			N/A	
VALUES IN PARENTHESIS ARE CALCULATED AND MAY CONTAIN ROUND OFF ERRORS			ITEM #	APPROX WEIGHT
			WFA1100	1.6 kg